WET WIPES MAXI AF

Stand: 161117

REF 100 766

Inhalt/Content/Contenu/Zawartość/Obsah/Tartalma:
150 Stück/pcs/unités/sztuk/ks/db

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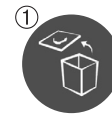
1 Processing of dispenser box
Remove lid from container

2 Wipe lid and inside of
container, the wipe tear-
off and reservoir with an
impregnated disinfection wipe
(observe exposure time of
disinfectant).

3 Allow container and
lid to dry

4 Cut the filled bag open at the
top and place in the container.
The disinfection wipes are
ready for use.

Type of application:

According to MPG: alcohol-based moist wipes for prophylactic surface disinfection of alcohol-resistant medical devices. Not suitable for the final (terminal) disinfection of invasive medical devices.

Notes on application:

Open the bag with wet wipes at the top and place in the intended dispensing tub. Pull the wipe from the middle of the fleece roll and thread from below through the star in the cover closure. Then press the cover firmly on the container and close tightly. Carefully wipe the areas to be disinfected and allow the disinfectant to work. Only to be used on alcohol-compatible materials (not suitable for acrylic glass). After removing a wipe, close the container carefully and tightly. After opening the bag, the wipes can only be used for 28 days in the corresponding dispensing container. The stick & lift label is removed from the bag, stuck to the dispensing tub and the use-by date added by hand. The dispensing tub is to be prepared according to the preparation instructions opposite with a suitable disinfection wipe. Wear suitable gloves during application.

Spectrum of activity:

Bactericidal (incl. MRSA) levurocidal (Candida albicans); tuberculocidal (M. terrae), "limited virucidal" (effective against enveloped viruses, such as, for example, HBV, HIV, HCV), and effective against adenoviruses, noroviruses and rotaviruses.

Exposure times (in relation to the impregnating liquid):

Bacteria according to EN 13727 and EN 13697 (high load):	Conc. 1 min
Candida albicans according to EN 13624 and EN 13697 (high load):	Conc. 1 min
TbB (tuberculosis) inactivation (M. terrae) acc. to EN 14348 (high load):	Conc. 1 min
Enveloped viruses (i.e. HBV, HIV, HCV) acc. to DVV (2012) with high load:	Conc. 1 min
Adenoviruses (EN 14476, high load):	Conc. 1 min
Noroviruses (EN 14476, high load):	Conc. 1 min
Rotaviruses (EN 14476, high load):	Conc. 1 min

Exposure times (in relation to the wet wipes):

Efficacy of the wet wipes acc. to EN 16615 with high load: Conc. 2 min
Bacteria and Candida albicans (tested according to the requirements of the VH, status 02.04.2015, high load): Conc. 2 min

Bag opened on: _____

Use-by date: _____ (= date of opening plus 28 days)

Signature: _____

Bacteria and Candida albicans (wipe impregnation system, included in VAH / DGHM listing):
Conc. 5 minutes (shortest exposure time for surface disinfectants in VAH list)

Contents: 150 pieces made of PES/VIS material mix; 20 cm x 30 cm impregnated.

Safety notes: H 228 inflammable solid. Keep away from heat/sparks/open flames/hot surfaces. No smoking. Wear protective gloves / eye protection. Store in a closed container.

100 g of impregnating liquid contain: 43.0 g ethanol, 8.0 g 1-propanol, 0.055 g didecyl-dimethyl-ammonium chloride.

Flash point acc. to DIN 51755: approx. 25°C. Application max. 50 ml per m². The total amount applied per room must not exceed 100 ml per m² of floor area. Keep container tightly sealed.

Only pass completely emptied containers to recycling. Dispose of residual material acc. to EAK 07 06 99.

Directive 648/2004/EC (impregnating liquid): Contains: < 5% cationic surfactants, fragrances

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